

# 2025

## Pressure Catalog

Enhancing quality  
of life is our  
driving force



**drive**

 **DeVilbiss®**  
HEALTHCARE



# Supportive Pressure Relief for Safety & Comfort

When it comes to products and solutions for individuals whose limited mobility puts them at risk of pressure injuries, we know that safety, comfort, and dependability are paramount. Drive DeVilbiss Healthcare has worked with clinical experts to develop effective products to meet the specialized needs and lifestyles of individuals who require dedicated support and care to minimize their risk of pressure injuries. Explore our array of pressure-reducing products, including cushions, pads, and a range of alternating pressure mattresses, and feel secure in the knowledge that you and your loved ones are in good hands with high-quality products from Drive.



# Pressure

## Premium Guard Gel Overlay

| Item # | Product Description    |
|--------|------------------------|
| 14893  | 34" W x 76" L x 3.5" H |
| 14901  | 42" W x 76" L x 3.5" H |

- Best in class gel overlay optimizes pressure redistribution and is now available in up to 500 lbs weight capacity
- 15 gel bladders cover entire patient contact area for the prevention and treatment of pressure ulcers
- Full length "head-to-toe" gel bladders minimize "bottoming out" and "pitting"
- New convoluted foam topper vastly improves patient comfort and air circulation
- 4 elastic corner straps secure overlay to mattress
- Urethane coated nylon cover is water resistant and vapor permeable
- Vinyl base is stain resistant and waterproof
- Bladders contain an aqueous non-toxic gel and polystyrene beads to minimize gel migration
- Meets flammability standard TB117

[Shop Now](#)



## Med-Aire Alternating Pressure Pump and Pad System

| Item #   | Product Description                |
|----------|------------------------------------|
| 14001E   | 36" W x 78" L - Variable Pressure  |
| 14002E   | 36" W x 78" L - Fixed Pressure     |
| 14001EFD | 36" W x 123" L - Variable Pressure |

- Quiet pump alternately inflates and deflates 130 individual bubble air cells
- 4 LPM (liters per minute) pump produces consistent air flow and pressure over a 5 minute cycle time and is available fixed or variable
- Variable pressure setting on pump allows for individualized comfort increasing patient compliance
- Durable, heavy gauge, premium vinyl bubble pad provides superior therapeutic treatment and is available with or without end flaps
- Pump has built-in brackets for mounting easily to bed

[Shop Now](#)



## Med-Aire Mattress Overlay Replacement

| Item # | Product Description |
|--------|---------------------|
|--------|---------------------|

|          |                |
|----------|----------------|
| 14003    | Pad Only       |
| 14003-EF | With End Flaps |

- 36" W x 78" L
- 36" W x 123" L

[Shop Now](#)

## Med-Aire Melody Alternating Pressure and Low Air Loss Mattress Replacement System

| Item # | Product Description |
|--------|---------------------|
|--------|---------------------|

|         |                           |
|---------|---------------------------|
| 14026   | Visual Alarm              |
| 14026-A | Visual Plus Audible Alarm |

- Designed to prevent, treat and heal pressure ulcers in the home or long term care setting
- Control unit offers 10 minute cycles when alternating and can also be used in static mode
- Adjustable patient weight settings allow for optimal immersion and patient comfort & compliance
- System includes 20 individual air cells to maximize immersion and envelopment with pillow feature for head stability
- 8" deep air cells are easily removed for convenient servicing
- Removable nylon cover is fluid resistant, low shear and vapor permeable
- CPR valve allows for rapid deflation
- Tubing quickly and easily disconnects

[Shop Now](#)

## Med-Aire 5" Alternating Pressure and Low Air Loss Overlay System

| Item #  | Product Description  |
|---------|----------------------|
| 14025NS | 36" W x 80" L x 8" H |

- Designed to prevent and treat pressure injuries in the home or long-term care setting, the Med-Aire 5" is an active mattress overlay system featuring 17 individual, 5" deep air cells. The 14025NS can be used on top of any standard foam or innerspring mattress and offers both alternating pressure and low air loss therapy to optimize pressure redistribution, shear and friction reduction, and microclimate control. The easy to use and intuitive control unit also offers a static feature and adjustable patient weight settings for optimal immersion, comfort, and compliance.
- Active mattress overlay system provides both alternating pressure and low air loss to optimize pressure redistribution, shear/friction reduction, and microclimate control
- Designed to prevent, treat and heal pressure injuries in the home or long term care setting
- Control unit offers 10 minute cycles when alternating and can also be used in static mode
- Adjustable patient weight settings allow for optimal immersion and patient comfort & compliance
- Visual alarm alerts when air pressure is low
- System includes 17 individual air cells to maximize immersion and envelopment
- Pillow feature maintains the 3 air cells at the head of the bed in static mode for head stability & patient comfort
- 5" deep air cells are easily removed for convenient servicing
- Removable Nylon cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- CPR valve allows for rapid deflation
- Tubing quickly and easily disconnects with EZ Lock/Quick Release design



[Shop Now](#)



## Med-Aire 8" Alternating Pressure and Low Air Loss Mattress System

| Item # | Product Description  |
|--------|----------------------|
| 14027  | 36" W x 80" L x 8" H |

- Active mattress replacement system provides both alternating pressure and low air loss to optimize pressure redistribution, shear/friction reduction, and microclimate control
- Designed to prevent, treat and heal pressure ulcers in the home or long term care setting
- Control unit offers 10 minute cycles when alternating and can also be used in static mode
- Adjustable patient weight settings allow for optimal immersion and patient comfort & compliance
- Audio and Visual alarm alerts when air pressure is low
- System includes 20 individual air cells to maximize immersion and envelopment with pillow feature for head stability & comfort
- 8" deep air cells are easily removed and replaced with air leakage plugs
- Removable Nylon stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- CPR valve allows for rapid deflation
- Tubing quickly and easily disconnects with EZ Lock/Quick Release design



[Shop Now](#)

## Med-Aire Plus 8" Alternating Pressure and Low Air Loss Mattress System

| Item #   | Product Description  |
|----------|----------------------|
| 14029    | 36" W x 80" L x 8" H |
| 14029-84 | 36" W x 84" L x 8" H |

- Active combination therapy system provides both alternating pressure and low air loss to optimize pressure redistribution, shear/friction reduction, and microclimate control
- Designed to prevent, treat and heal pressure ulcers in the home or long term care setting
- Digital control unit offers 4 cycle times and adjustable comfort weight settings for optimal immersion and patient satisfaction
- 3 diagnostic audio and visual alarms alert for issues with air pressure, the system, or power
- System includes 20 individual air cells to maximize immersion and envelopment with pillow feature for head stability & comfort
- Advanced Clinical and Care Giver Features: Auto-Firm, Seat Inflate, microprocessing feedback, lock-out and more!
- Cell-on-cell designed air cells prevent "bottoming out" and provide up to 24 hours of power outage protection
- 8" deep air cells are easily removed with quick disconnects for convenient servicing
- Removable 4 way stretch PU cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- CPR valve allows for rapid deflation
- Tubing quickly and easily disconnects with EZ Lock/Quick Release design



[Shop Now](#)

## Med-Aire Plus 8" Alternating Pressure and Low Air Loss Mattress System with 10" Defined Perimeter

| Item #  | Product Description |
|---------|---------------------|
| 14029DP | 36"W x 80"L x 10"H  |

- Active combination therapy system provides both alternating pressure and low air loss to optimize pressure redistribution, shear/friction reduction, and microclimate control
- Designed to prevent, treat and heal pressure ulcers in the home or long term care setting
- Digital control unit offers 4 cycle times and adjustable comfort weight settings for optimal immersion and patient satisfaction
- 3 diagnostic audio and visual alarms alert for issues with air pressure, the system, or power
- 10" deep static perimeter tubing surrounds the length of the mattress preventing patient falls or entrapment
- Pillow feature provides head stability & comfort
- Advanced Clinical and Care Giver Features: Auto-Firm, Seat Inflate, microprocessing feedback, lock-out and more!
- Cell-on-cell designed air cells prevent "bottoming out" and provide up to 24 hours of power outage protection
- 8" deep air cells are easily removed with quick disconnects for convenient servicing
- Removable 4 way stretch PU cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- CPR valve allows for rapid deflation
- Tubing quickly and easily disconnects with EZ Lock/Quick Release design



[Shop Now](#)



## Med-Aire Plus 10" Bariatric Alternating Pressure and Low Air Loss Mattress Replacement System

| Item # | Product Description   |
|--------|-----------------------|
| 14030  | 42" W x 80" L x 10" H |
| 14048  | 48" W x 80" L x 10" H |
| 14054  | 54" W x 80" L x 10" H |
| 14060  | 60" W x 80" L x 10" H |

- Active combination bariatric therapy system provides both alternating pressure and low air loss to optimize pressure redistribution, shear/friction reduction, and microclimate control
- Designed to prevent, treat and heal all stages of pressure ulcers in the home or long term care setting
- Digital control unit offers 4 cycle times and adjustable comfort weight settings for optimal immersion and patient satisfaction
- 3 diagnostic audio and visual alarms alert for issues with air pressure, the system, or power
- System includes 20 individual air cells to maximize immersion and envelopment with pillow feature for head stability & comfort
- Advanced Clinical and Care Giver Features: Auto-Firm, Seat Inflate, auto-recovery, microprocessing feedback, lock-out and more!
- Cell-on-cell designed air cells prevent "bottoming out" and provide up to 24 hours of power outage protection
- 10" deep air cells are easily removed with quick disconnects for convenient servicing
- Removable quilted Nylon cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- CPR valve allows for rapid deflation
- Tubing quickly and easily disconnects with EZ Lock/Quick Release design



[Shop Now](#)

## Med-Aire Assure 5" Air + 3" Foam Base Alternating Pressure and Low Air Loss Mattress System

| Item # | Product Description |
|--------|---------------------|
| 14530  | 36"W x 80"L x 8"H   |

- Active mattress replacement system provides both alternating pressure and low air loss to optimize pressure redistribution, shear/friction reduction, and microclimate control
- Designed to prevent, treat and heal pressure ulcers in the home or long term care setting
- Control unit offers 10 minute cycles when alternating and can also be used in static mode
- Adjustable patient weight settings allow for optimal immersion and patient comfort & compliance
- Audio and Visual alarm alerts when air pressure is low
- System includes 17 individual air cells over a 3" foam base for power outage protection
- Pillow feature provides head stability & comfort
- 5" deep air cells are easily removed for convenient servicing
- Removable quilted Nylon cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- CPR valve allows for rapid deflation
- Tubing quickly and easily disconnects with EZ Lock/Quick Release design



[Shop Now](#)

## PreserveTech™ Med-Aire Edge Alternating Pressure & Low Air Loss Mattress Replacement System

| Item #  | Product Description            |
|---------|--------------------------------|
| 14360-B | 36" W x 80" L x 8" H - Analog  |
| 14360-P | 36" W x 80" L x 8" H - Digital |

- Active mattress replacement system provides both alternating pressure and low air loss to optimize pressure redistribution, shear/friction reduction, and micro-climate control
- Designed to help prevent, treat and heal pressure injuries
- Control unit offers 10-minute cycles when alternating and can also be used in static mode (exclusive to #14360-B)
- Control dial to optimize patient comfort levels by weight (exclusive to #14360-B)
- Advanced digital control unit offers 4 adjustable APM cycle times (10/15/20/25), and adjustable comfort weight settings for optimal immersion & patient satisfaction (exclusive to #14360-P)
- 9 deluxe clinical and care giver features: auto-firm, seat inflate, micro processing feedback, lock-out and more (exclusive to #14360-P)
- Audio and Visual alarm alerts when air pressure is low
- Digital control unit features 3 diagnostic audio and visual alarms to alert for issues with air pressure, the system, or power (exclusive to #14360-P)
- System includes 20 individual air cells to maximize immersion and envelopment and redistribute pressure with pillow feature for head stability and comfort.
- 8" firm foam perimeter for safe ingress and egress (getting into and out of the bed), patient transfers and power outage protection.
- 3" foam base for power outage protection
- 5" deep air cells are easily removed with quick disconnects for convenient servicing
- Removable and machine washable quilted 4-way stretch PreserveTech™ cover is fluid resistant, low shear and vapor permeable
- CPR valve allows for rapid deflation
- Tubing quickly and easily disconnects with EZ Lock/Quick Release design
- 45db decibel level



[Shop Now](#)

## PreserveTech™ Harmony "True" Low Air Loss Tri-Therapy System

| Item #    | Product Description   |
|-----------|-----------------------|
| 14200N    | 36" W x 80" L x 10" H |
| 14200N-42 | 42" W x 80" L x 10" H |

- Designed to prevent, treat and heal pressure injuries in the home or long term care setting.
- Advanced, user friendly high-flow digital control unit offers (4) therapy cycle times ranging from 5 to 20 minutes, and (10) adjustable comfort weight settings for optimal immersion and patient satisfaction.
- High air flow "true" low air loss system wicks away moisture and humidity providing optimal micro-climate control
- Features three diagnostic audible and visual alarms to alert for issues with air pressure, power or the system.
- 20 individual, 8" deep air cells provide optimal patient comfort and pressure redistribution.
- 2" foam base offers protection in the event of a power outage or system failure.
- The cell material, Nylon/TPU blend, provides a specialty surface that conforms to the specific shape of the patient (immersion), minimizing soft tissue distortion, reducing bone and muscle friction.
- Quick cell disconnects offer convenient servicing and auto firm feature inflates air cells in under two minutes for rapid set up and ease of care giving.
- Removable, quilted four way stretch cover is manufactured with an anti-microbial\* agent and is fluid resistant, low shear and vapor permeable helping protect the skin
- Built in cable management with (5) side mattress loops to keep cables off the floor, reducing tripping hazards.
- Snap fit buckle allows the system to be compactly wrapped for quick and easy storage within the included perfect fit carry bag
- CPR connector in case of emergency procedures releases the air immediately from the mattress for rapid deflation
- Tubing quickly and easily disconnects with EZ lock/quick release design
- Clinical and caregiver features include pulsation therapy, auto-firm, micro-processing feedback, lock-out and more!



[Shop Now](#)



## PreserveTech™ Lateral Rotation System With On Demand Low Air Loss

| Item #  | Product Description   |
|---------|-----------------------|
| LS9500N | 36" W x 80" L x 10" H |

- Active mattress replacement system incorporating the functions of lateral rotation and low air loss to assist with pulmonary blood flow and mucosal drainage aiding in the prevention and treatment of pressure injuries in a single portable system.
- Combination therapy active system incorporates the functions of lateral rotation and low air loss to assist with pulmonary blood flow and mucosal drainage, as well as provide microclimate control to keep the patient cool and dry aiding in the prevention and treatment of pressure injuries in a single complete system.
- Mattress features 20, 8" individual split air cells provide customizable turn and rotation angles for optimal patient comfort and pressure redistribution.
- 2" foam base offers protection in the event of a power outage or system failure.
- 10" high side air bolsters provide fall prevention protection and facilitate safe patient transfers.
- Unique split air cells provide customizable turn time and rotation angles for patient comfort and satisfaction.
- Control unit offers 10 comfort levels, 10,20, 30 & 60 minute adjustable turning cycles and maximum turn angle adjustments.
- Turning modes include combinations of ¼, ½, ¾, and adjustable turning positions up to 40 degrees with left, right, both or none (static) directions.
- Removable quilted four way stretch polyurethane coated cover is manufactured with an anti-microbial\* agent and is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture.
- Pillow feature provides head stability and comfort.
- CPR connector allows for rapid deflation.
- Tubing quickly and easily disconnects with EZ lock/quick release design.
- Built in side loop cable management system reduces tripping hazards.
- Audio and Visual alarm alerts when air pressure is low.
- Advanced Clinical and Care Giver Features: Auto-Firm, auto-recovery, microprocessing feedback, lock-out and more!



[Shop Now](#)

## Balanced Aire Powered Self Adjusting Convertible Mattress

| Item #       | Product Description                |
|--------------|------------------------------------|
| BA9600-NP    | 35" W x 80" L x 7" H - Non-Powered |
| BA9600-NP-84 | 35" W x 84" L x 7" H - Non-Powered |
| BA9600-P     | 36" W x 80" L x 7" H - Powered     |
| BA9600-P-84  | 36" W x 84" L x 7" H - Powered     |
| BA9600-NP-42 | 42" W x 80" L x 7" H - Non-Powered |
| BA9600-P-42  | 42" W x 80" L x 7" H - Powered     |

- Versatile reactive mattress replacement system with self adjusting air/foam cells provides outstanding pressure redistribution and shear/friction reduction
- Designed to prevent, treat and heal pressure ulcers in the home or long term care setting
- Add an alternating pressure pump and the Balanced Aire easily converts to a dynamic powered system (14005E or BA9600-PUMP)
- 2.5" Hybrid Visco foam top layer provides immersion, premium comfort and heat dissipation
- High density foam perimeter provides firm edge support and facilitates safe patient transfers
- Heel slope redistributes pressure away from vulnerable heels
- Self-adjustment can be open cell or closed cell to meet individual needs
- Self-Adjusting air cells recognize body weight and pressure and adjust to reduce interface pressure and improve patient comfort
- Eight non-circular cylinders prevent "rolling" or "gapping"
- Removable 4 way stretch PU cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Arrives patient ready easy to deliver, set up and transport
- Foam and cover are antimicrobial, antifungal, and latex free
- Meets flammability standards CFR1632 and CRF1633



[Shop Now](#)

## Quick Fill Pump With Adapters

**Item #**

14427QF

- The Quick Fill Pump is designed for fast inflation of pressure prevention air mattresses, ensuring a rapid and efficient setup. It includes four adapters for compatibility with most pressure prevention mattress replacement systems.



[Shop Now](#)

## Therapeutic 5 Zone Support Mattress

**Item #**

15019

**Product Description**

36" W x 80" L

- High quality latex free therapeutic foam mattress provides pressure redistribution and shear/friction reduction
- Deluxe horizontal, cross-cut foam mattress provides comfort, support, and pressure redistribution over 5 therapeutic pressure zones
- Non skid bottom keeps mattress safely in place
- Removable Nylon cover with welded seams is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Newly added handles make setup, delivery, and cleaning easier
- Intuitive tie-down straps easily fasten to virtually any bed frame for a secure fit to reduce mattress migration and falls
- Meets flammability standard CFR 1632 and CFR 1633




[Shop Now](#)

## Gravity 6 Long Term Care Pressure Redistribution Mattress

| Item # | Product Description |
|--------|---------------------|
| 15670  | 36" W x 80" L x 6"  |

- High quality latex free therapeutic foam mattress provides pressure redistribution and shear/friction reduction
- Intuitive tie-down straps easily fasten to virtually any bed frame for a secure fit to reduce mattress migration and falls
- Full length articulation cuts to extend mattress life cycle and improve conformity to all bed positions
- Newly added handles make setup, delivery, and cleaning easier
- Deluxe horizontal, cross-cut foam mattress provides comfort, support, and pressure redistribution over 5 therapeutic pressure zones
- Better compliance and durability than innerspring mattresses allows for improved patient satisfaction
- Non-skid vinyl bottom cover keeps mattress safely in place
- Heel section is sloped offloading pressure
- Removable multi-stretch cover with welded seams is fluid resistant, low shear and vapor permeable helping protect the skin
- Softer and longer-lasting durable cover tolerates disinfectant cleaning supporting infection control, refer to user manual for suggested disinfectants
- Meets flammability standard CFR 1632 and CFR 1633

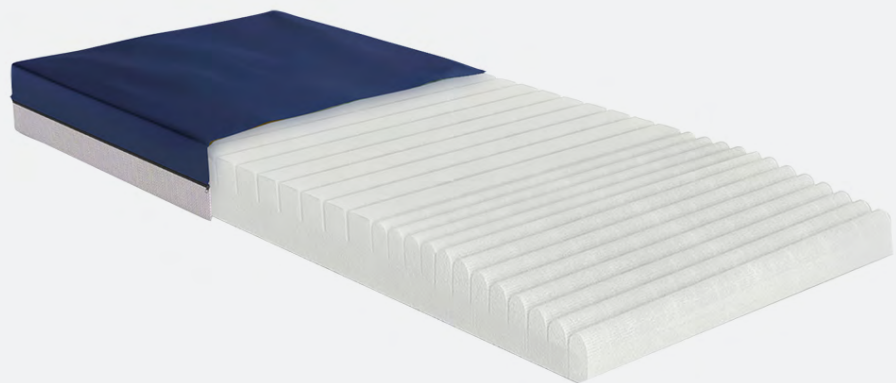


[Shop Now](#)



## Gravity 7 Long Term Care Pressure Redistribution Mattress

| Item # | Product Description                         |
|--------|---|
| 15770  | 36" W x 80" L x 6" H - Without Raised Rails |
| 15777  | 36" W x 80" L x 6" H - With Raised Rails    |
| 15784  | 36" W x 84" L x 6" H - Without Raised Rails |
| 15785  | 36" W x 84" L x 6" H - With Raised Rails    |
| 15876  | 36" W x 76" L x 6" H - Without Raised Rails |
| 15877  | 36" W x 76" L x 6" H - With Raised Rails    |

- Therapeutic die-cut foam spans the full width for maximum patient contact
- Drive's most popular therapeutic foam mattress maximizes value and comfort with 3 full layers of latex free foam
- Full length articulation cuts to extend mattress life cycle and improve conformity to all bed positions
- Bottom layer provides full length horizontally scored articulation cuts extending the durability and life of the mattress
- High-density foam perimeters remain beneath for safe patient transfers
- Available with or without raised side rails
- Visco foam heel section offloads pressure from vulnerable heels
- Removable multi-stretch cover with welded seems is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Softer, more supple, longer lasting cover fabric for enhanced comfort
- Reduce mattress migration and reduce falls- easy to use and intuitive "d" ring straps quickly and easily secure mattress to virtually any bed frame
- High-density foam perimeters remain beneath for safe patient transfers
- New handles make setup, delivery, and cleaning easier
- Unparalleled non-prorated 7-year warranty with Date Stamped Foam
- Meets flammability standard CFR 1632 and CFR 1633

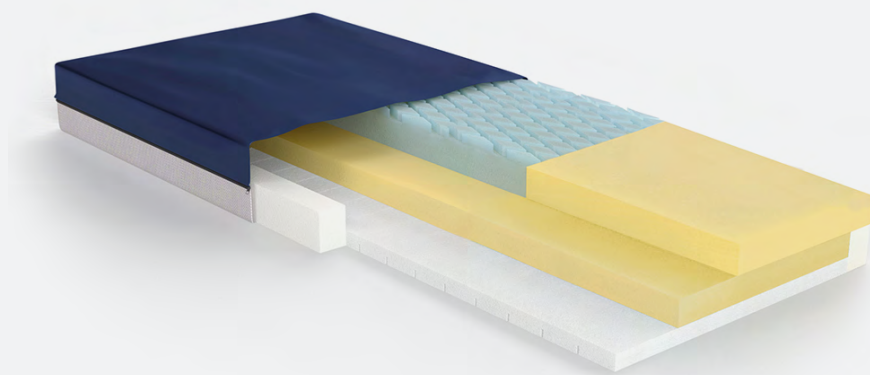


[Shop Now](#)

## Gravity 8 Deluxe Long Term Care Pressure Redistribution Mattress

| Item # | Product Description                         |
|--------|---|
| 15870  | 36" W x 80" L x 6" H - Without Raised Rails |
| 15884  | 36" W x 84" L x 6" H - Without Raised Rails |
| 15886  | 36" W x 76" L x 6" H - Without Raised Rails |
| 15888  | 36" W x 80" L x 6" H - With Raised Rails    |

- Therapeutic die-cut foam spans the full width for maximum patient contact
- High end therapeutic foam mattress maximizes value and comfort with 3 full layers of latex free foam
- Full length articulation cuts to extend mattress life cycle and improve conformity to all bed positions
- Bottom layer provides full length horizontally scored articulation cuts extending the durability and life of the mattress
- High-density foam perimeters remain beneath for safe patient transfers
- Available with or without raised side rails
- Newly added handles make setup, delivery, and cleaning easier
- Heel slope with channel cut foam reduces pressure on patient heels
- Removable multi-stretch cover with welded seems is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Meets flammability standard CFR 1632 and CFR 1633

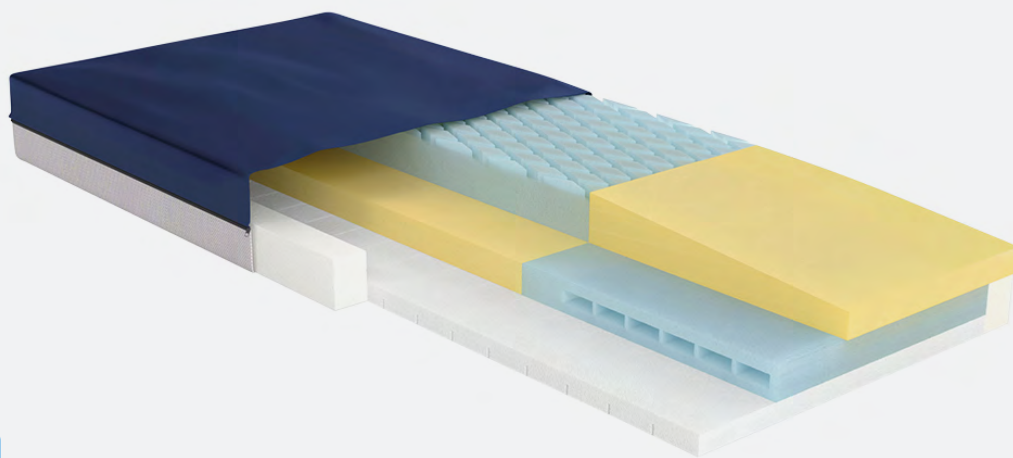


[Shop Now](#)

## Gravity 9 Premium Long Term Care Pressure Redistribution Mattress

| Item #     | Product Description                         |
|------------|---|
| 15970      | 36" W x 80" L x 6" H - Without Raised Rails |
| 15977      | 36" W x 80" L x 6" H - With Raised Rails    |
| 15984      | 36" W x 84" L x 6" H - Without Raised Rails |
| 15996      | 36" W x 76" L x 6" H - Without Raised Rails |
| 15970-4280 | 42" W x 80" L x 6" H - Without Raised Rails |
| 15970-4284 | 42" W x 84" L x 6" H - Without Raised Rails |
| G9-39X80X6 | 39" W x 80" L x 6" H - Multi-Zone           |

- Therapeutic die-cut foam spans the full width for maximum patient contact
- Head to toe hybrid visco memory foam top layer
- Drive's premium therapeutic Gravity foam mattress maximizes value and comfort with multiple layers of latex free foam
- Full length articulation cuts to extend mattress life cycle and improve conformity to all bed positions
- Full top layer of open cell, visco-elastic memory foam improves immersion and pressure redistribution
- Bottom layer provides full length horizontally scored articulation cuts extending the durability and life of the mattress
- High-density foam perimeters remain beneath for safe patient transfers
- Available with or without raised side rails
- Newly added handles make setup, delivery, and cleaning easier
- Heel slope with channel cut foam reduces pressure on patient heels
- Removable multi-stretch cover with welded seems is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Meets flammability standard CFR 1632 and CFR 1633

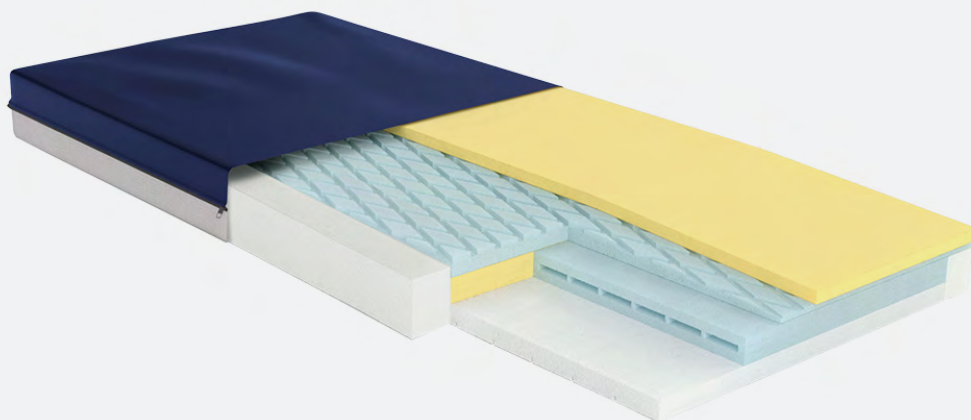


[Shop Now](#)

## Multi-Ply 6500 Dynamic Elite Pressure Redistribution Foam Mattress

| Item #       | Product Description  |
|--------------|----------------------|
| 6500-DE-2-FB | 36" W x 80" L x 6" H |

- Therapeutic die-cut foam spans the full width for maximum patient contact
- Full length articulation cuts to extend mattress life cycle and improve conformity to all bed positions
- High-density foam perimeters remain beneath for safe patient transfers
- Available with or without raised side rails
- Heel section is sloped with patented LuraQuilt foam offloading pressure from vulnerable heels
- Newly added handles make setup, delivery, and cleaning easier
- Removable premium, multi-stretch cover with welded seams is fluid resistant, low shear and vapor permeable helping protect the skin
- Meets flammability standard CFR 1632 and CFR 1633

[Shop Now](#)

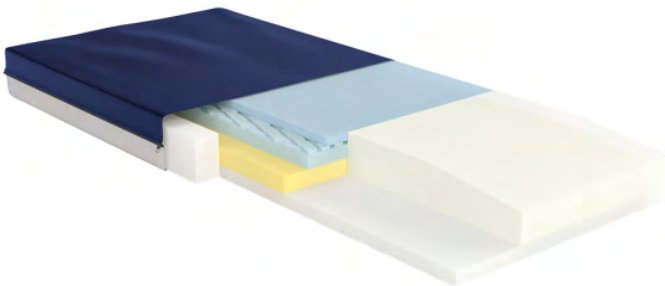

## Multi-Ply ShearCare 1500 Pressure Redistribution Foam Mattress

| Item #       | Product Description |
|--------------|---------------------|
| 1500SC-42-FB | 42" W x 80" L       |
| 1500SC-48-FB | 48" W x 80" L       |

- Bariatric, latex free, therapeutic foam mattress provides pressure redistribution and shear/friction reduction
- Full width top layer is high-resilient, high-density die-cut foam for pressure redistribution and maximum support
- Intuitive tie-down straps easily fasten to virtually any bed frame for a secure fit to reduce mattress migration and falls
- Premium foam heel section offloads pressure from vulnerable heels
- Breathable stretch cover with welded seams for better water resistance
- Newly added handles make setup, delivery, and cleaning easier
- Meets flammability standard CFR 1632 and CFR 1633

[Shop Now](#)






## Quick 'N Easy Comfort Mattress

| Item # | Product Description |
|--------|---------------------|
| 15076  | 36" W x 80" L x 6"  |

- Durable high quality latex free polyurethane foam mattress for utmost comfort and support
- Upgraded handles make setup, delivery, and cleaning easier
- Full length articulation cuts to extend mattress life cycle and improve conformity to all bed positions
- Intuitive tie-down straps easily fasten to virtually any bed frame for a secure fit to reduce mattress migration
- Compression technology creates a 14" x 14" x 37.5" box footprint creating convenience for transport and delivery
- Softer and longer-lasting removable cover with welded seams is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Meets flammability standard CFR 1632 and CFR 1633

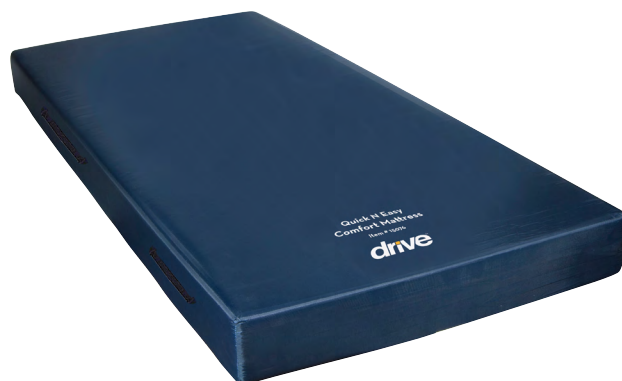[Shop Now](#)

## Defined Perimeter Mattress Cover

| Item #   | Product Description          |
|----------|------------------------------|
| 14333    | 36"-39" W x 80" L x 6"-8" H  |
| 14333-42 | 39"-43" W x 80" L x 6"-10" H |

- The Defined Perimeter Mattress Cover creates a raised rail, defined perimeter for enhanced fall prevention, without using patient restraints
- Open in mid-section for safe and easy patient egress and ingress
- Built-in medium density foam bolsters safely and comfortably define edges of bed for patient and prevent falling
- Easy-to-install and secure
- Breathable Nylon makes this a truly universal cover, compatible with both foam and powered air surfaces

[Shop Now](#)

## 3" Foam Retail Cushion

**Item #**  
RTL14910

- 3" foam cushion in attractive retail package with cutout to feel cover
- Comes with durable, machine washable, jacquard zippered cover
- Weight Capacity 250 lbs
- Made of soft, premium quality foam
- Versatile use for any location where seated comfort is needed



[Shop Now](#)

## Premier One General Use Foam Cushion

| Item #   | Product Description  |
|----------|----------------------|
| M8061-PB | 14" W x 14" L x 3" H |
| M8062-PB | 14" W x 16" L x 3" H |
| M8063    | 16" W x 16" L x 3" H |
| M8064    | 18" W x 16" L x 3" H |
| M8065    | 20" W x 16" L x 3" H |
| M8070    | 20" W x 18" L x 3" H |
| M8077    | 22" W x 18" L x 3" H |

- The Premier One General Use Foam cushion is designed to provide comfort and stability with a simple, economical design
- Aerated foam helps keep the user cool and dry
- The low-shear water-resistant stretch cover helps maintain sensitive skin integrity
- Convenient 12" tie straps secure the cushion to the wheelchair
- The waterproof vinyl bottom protects the wheelchair seat
- Meets flammability standard TB117

[Shop Now](#)



## Molded Foam General Use Wheelchair Cushion

| Item # | Product Description  |
|--------|----------------------|
| 14880  | 16" W x 16" L x 2" H |
| 14887  | 18" W x 16" L x 2" H |
| 14881  | 20" W x 16" L x 2" H |
| 14908  | 18" W x 18" L x 2" H |
| 14909  | 20" W x 18" L x 2" H |

- A unique construction of high-resilient molded foam provides optimal support, stabilization and pressure redistribution
- Designed for low risk patients seeking comfort and support
- Removable, zippered, low-shear stretch cover is vapor permeable and water-resistant
- Cover has a "Non-Slip" vinyl bottom for comfort and safety
- Convenient 12" tie straps secure the cushion to wheelchair
- Meets flammability standard TB117

[Shop Now](#)

## Lumbar Support General Use Wheelchair Back Cushion

| Item # | Product Description    |
|--------|------------------------|
| 14889  | 18" W x 17" L x 2.5" H |
| 14906  | 16" W x 17" L x 2.5" H |
| 14920  | 20" W x 17" L x 2.5" H |
| 8033   | 22" W x 19" L x 2.5" H |

- Memory-foam back cushion provides support, comfort and a tactile curve to improve posture
- Durable Masonite board prevents hammocking in sling-back wheelchair Easily removed as per patient preference
- Conforms to lumbar curve to support lower back
- Sealed liner protects foam core from moisture
- Nylon stretch cover is water resistant, non-slip vinyl back prevents the cushion from sliding and is waterproof
- Adjustable Straps easily attach to wheelchair
- Meets flammability standard TB117

[Shop Now](#)

## Adjustable Tension General Use Wheelchair Back Cushion

| Item # | Product Description |
|--------|---------------------|
| 14300  | 16"-21" Adjustable  |
| 14301  | 22"-26" Adjustable  |

- Versatile back support system that improves patient comfort and posture
- 1" thick high density foam pad provides optimal pressure redistribution
- Secures directly over wheelchair canes for added strength, stability and safety
- 5 Adjustment straps featured for multi-purpose back support and customization to individual seating needs
- Handy back pocket is convenient for carrying keys, water bottle, medications, tissues and more
- Breathable cover is fluid resistant and easy to clean
- Cushion folds with the wheelchair for easy transport and storage
- Meets flammability standard TB117

[Shop Now](#)

## Wedge Cushion with Stretch Cover

| Item # | Product Description     |
|--------|-------------------------|
| 8134   | 18" W x 16" D x 4"-2" H |

- Wedge design helps prevent patient from sliding forward
- Firm foundation reduces pressure caused by sling-seat hammocks
- The top cover is fluid resistant, low shear stretch fabric
- The bottom cover is waterproof
- High-density foam increases stability and comfort

[Shop Now](#)





## General Use Wheelchair Cushion Kits

| Item #   | Product Description       |
|----------|---------------------------|
| 14880KIT | Fits 16" Wide Wheelchairs |
| 14887KIT | Fits 18" Wide Wheelchairs |
| 14881KIT | Fits 20" Wide Wheelchairs |

- Kit includes one seat and one back cushion each for added comfort & support
- Seat cushion is Drive's Molded Foam General Use cushion for optimal support and stabilization
- Back cushion is Drive's Lumbar Support General Use cushion for supporting the lower back and improving posture
- Meet flammability standard TB117

[Shop Now](#)



## Gel "E" 2" General Use Gel/Foam Wheelchair Cushion

| Item # | Product Description  |
|--------|----------------------|
| 14888  | 18" W x 16" L x 2" H |

- The Medicare approved General Use 2" gel-foam wheelchair cushion will provide comfort and support for the wheelchair bound patient
- High-Density foam encapsulates a viscous gel bladder that provides proper immersion and envelopment for pressure redistribution
- The cover is fluid resistant low-shear nylon zippered and removable with a durable vinyl bottom
- Carry handle included for convenient transport
- Meets flammability standard TB117

[Shop Now](#)



## Gel-U-Seat™ Lite General Use 2" Gel/Foam Wheelchair Cushion

| Item # | Product Description  |
|--------|----------------------|
| 8040-1 | 16" W x 16" L x 2" H |
| 8040-2 | 18" W x 16" L x 2" H |
| 8040-3 | 20" W x 16" L x 2" H |
| 8040-5 | 18" W x 18" L x 2" H |
| 8040-4 | 20" W x 18" L x 2" H |
| 8040-7 | 22" W x 18" L x 2" H |
| 8040-9 | 24" W x 18" L x 2" H |

- Offering great versatility, the Gel-U-Seat™ Lite provides comfort to almost all wheelchair users
- Unique dual-chambered gel bladder keeps the pelvis in contact with pressure redistributing cooling gel
- The cover is fluid resistant low-shear nylon zippered and removable with a durable vinyl bottom
- Secures easily to wheelchair with 12" tie straps
- Meets flammability standard TB117

[Shop Now](#)



## Gel-U-Seat™ Skin Protection 3" Gel/Foam Wheelchair Cushion

| Item # | Product Description  |
|--------|----------------------|
| 8046-1 | 16" W x 16" L x 3" H |
| 8046-2 | 18" W x 16" L x 3" H |
| 8046-3 | 20" W x 16" L x 3" H |
| 8046-4 | 20" W x 18" L x 3" H |
| 8046-5 | 22" W x 18" L x 3" H |
| 8046-7 | 24" W x 18" L x 3" H |

- The Gel-U-Seat™ skin protection cushion has a high density foam shell surrounding a uniquely designed gel bladder
- Dual chambered gel bladder provides proper immersion and envelopment for pressure redistribution
- Cover is water resistant and low shear stretch to reduce heat and moisture build-up
- Convenient 12" tie straps secure the cushion to the wheelchair
- Meets flammability standard TB117

[Shop Now](#)



## Gel "E" Skin Protection 3" Gel/Foam Wheelchair Cushion

| Item # | Product Description  |
|--------|----------------------|
| 14891  | 16" W x 16" L x 2" H |
| 14886  | 18" W x 16" L x 3" H |
| 14878  | 18" W x 18" L x 3" H |
| 14894  | 20" W x 16" L x 3" H |
| 14902  | 22" W x 18" L x 3" H |


- Pressure redistributing cushion designed to assist in the prevention, treatment and management of pressure injuries
- Viscous gel bladder is encased in polyurethane foam shell and provides optimal pressure redistribution, support and comfort
- Cover top is a urethane-coated nylon that is low-shear, vapor permeable and water-resistant
- Cover base is waterproof vinyl for durability
- Carry handle included for convenient transport
- Meets flammability standard TB117

[Shop Now](#)

## Titanium Skin Protection and Positioning Gel/Foam Wheelchair Cushion

| Item #   | Product Description    |
|----------|------------------------|
| FPT-1    | 16" W x 16" L x 3.5" H |
| FPT-2    | 18" W x 16" L x 3.5" H |
| FPT-1818 | 18" W x 18" L x 3.5" H |
| FPT-2018 | 20" W x 18" L x 3.5" H |
| FPT-4    | 22" W x 18" L x 3.5" H |

- The Titanium cushion is specifically designed to offer superior pressure redistribution and positioning
- High-resilient, high-density foam combined with a dual chamber gel bladder encourages proper patient posture through immersion
- Middle contour acts as a built-in leg abductor
- Low-shear, fluid-resistant stretch cover helps to maintain sensitive skin tissue
- Non-skid bottom offers additional stability
- Meets flammability standard TB117

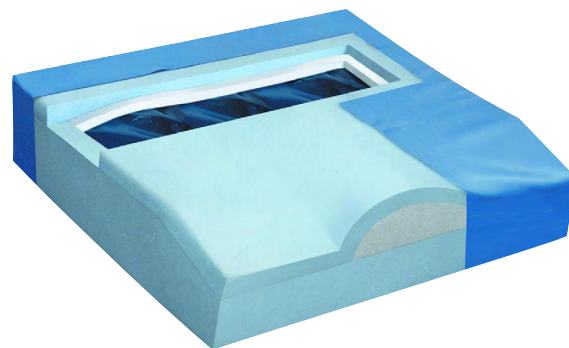[Shop Now](#)

## Balanced Aire Adjustable Skin Protection Cushion

| Item #      | Product Description  |
|-------------|----------------------|
| 8047-16     | 16" W x 16" L x 4" H |
| 8047-16-2   | 16" W x 16" L x 2" H |
| 8047-18     | 18" W x 16" L x 4" H |
| 8047-1818   | 18" W x 18" L x 4" H |
| 8047-1818-2 | 18" W x 18" L x 2" H |
| 8047-18-2   | 18" W x 16" L x 2" H |
| 8047-20     | 20" W x 16" L x 4" H |
| 8047-2018   | 20" W x 18" L x 4" H |
| 8047-20-2   | 20" W x 16" L x 2" H |
| 8047-2218   | 22" W x 18" L x 4" H |

- Flexible interconnected, independent air cells allow for deep immersion into the cushion
- Ideal for clients with limited seated tolerance, or those at risk for pressure injuries
- Low-shear stretch cover with adjustable locking straps and non-skid bottom is easily removed and washed
- Hand pump included to adjust the cushion
- Repair kit included for quick fixes
- Meets flammability standard TB117

[Shop Now](#)



## Bed Wedges (Cloth Cover)

| Item # | Product Description   |
|--------|-----------------------|
| 3826   | 23" W x 23" L x 10" H |
| 3827   | 23" W x 23" L x 10" H |

- Bed wedges with a cloth cover are designed to help elevate and support the neck
- Helps to relieve discomfort associated with respiratory problems, poor circulation and back or neck pain
- 100% polyurethane foam

[Shop Now](#)



## Folding Bed Wedges

| Item #  | Product Description   |
|---------|-----------------------|
| RTL3825 | 23" W x 23" L x 7" H  |
| RTL3826 | 23" W x 23" L x 10" H |
| RTL3827 | 23" W x 23" L x 12" H |

- Designed to elevate and support the neck as well as help alleviate respiratory conditions
- Helps to relieve discomfort from back or neck pain and assists with poor circulation
- Helps to alleviate discomfort from pregnancy
- Great for sleeping with acid reflux or heartburn
- Compact folding design allows for easy storage and the carry handle is perfect for travel
- Available in 3 sizes

[Shop Now](#)

## Compressed Bed Wedge Cushion

**Item #**  
RTL1490COM

- Designed to elevate and support the neck as well as help alleviate respiratory conditions
- Helps to relieve discomfort from back or neck pain and assists with poor circulation
- Great for sleeping with acid reflux or heartburn
- Cushion automatically expands to full size when removed from packaging
- Comfortable and durable foam construction
- Removable, machine-washable cover

[Shop Now](#)



## Foam Ring Cushion

**Item #**

RTLPC23388

**Product Description**

16" W x 12.75" L x 2.75" H

- Provides excellent comfort and support while seated
- Open center is designed to provide maximum pressure relief
- Comes with durable, washable, zippered jacquard cover
- Made from soft, premium quality therapeutic foam
- Weight Capacity: 250lbs

[Shop Now](#)

## Compressed Coccyx Cushion

**Item #**

RTL1491COM

- Cut-out reduces pressure on spine and coccyx (tail bone) when seated
- Cushion's tilt restores spines' natural "S" shape
- Great for home, office and travel
- Cushion automatically expands to full size when removed from packaging
- Comfortable and durable foam construction
- Removable, machine-washable cover

[Shop Now](#)

## Compressed Foam Ring

| Item #     | Product Description        |
|------------|----------------------------|
| RTL1492COM | 16.5" W x 12.5" L x 3.5" H |

- Specially designed to reduce pressure on sensitive areas when sitting for extended periods
- Unique shape offers relief of coccyx (tail bone) pain, hemorrhoids and post-surgery discomfort
- Donut shape comfortably conforms to body contours
- Cushion automatically expands to full size when removed from packaging
- Comfortable and durable foam construction
- Removable, machine-washable cover

[Shop Now](#)

## Compressed Lumbar Cushion

| Item #     | Product Description    |
|------------|------------------------|
| RTL1493COM | 13.5" W x 13" L x 4" H |

- Transforms an ordinary chair by providing ergonomic comfort and back support
- Ideal for those who spend extended periods sitting
- Subtle contour promotes proper posture while offering support and minimizing back pain
- Attached strap helps keep cushion securely in place
- Great for home, office and travel—even in the car
- Cushion automatically expands to full size when removed from packaging
- Comfortable and durable foam construction
- Removable, machine-washable cover

[Shop Now](#)

## Compressed Posture Support Cushion

| Item #     | Product Description     |
|------------|-------------------------|
| RTL1494COM | 16.4" W x 8" L x 4.5" D |

- Transforms an ordinary chair by providing ergonomic comfort and back support
- Firm yet plush for superior posture support
- Helps reduce discomfort in the head, neck, shoulders and back
- Ergonomically designed to promote proper spinal alignment
- Attached strap helps keep cushion securely in place
- Great for home, office and travel—even in the car
- Also great for use behind the knees and ankles
- Cushion automatically expands to full size when removed from packaging
- Comfortable and durable foam construction
- Removable, machine-washable cover

[Shop Now](#)

## Comfort Touch™ Knee Support Cushion

| Item #     | Product Description   |
|------------|-----------------------|
| RTL2017CTK | 10" W x 7.5" L x 8" H |

- Alleviates discomfort by promoting proper alignment of spine, knees and hips
- Contoured shape fits comfortably between the legs with adjustable, secure strap to keep cushion in place
- Soft velour cover enhances relaxation and is easily removable for cleaning
- Ideal for side sleepers, people with lower back or Sciatica nerve pain, or women who are pregnant

[Shop Now](#)

## Comfort Touch™ Heated Lumbar Support

| Item #     | Product Description        |
|------------|----------------------------|
| RTL2017CTL | 12.6" W x 13.2" L x 4.1" H |

- Single heat setting provides continuous heat therapy to ease back pain and muscle tension
- Ergonomically designed to provide postural support
- Soft velour cover provides luxurious feel for added comfort and is easily removable for cleaning
- Can be used on your seat of choice wherever a power supply is available

[Shop Now](#)



## Comfort Touch™ Neck Support Pillow

| Item #     | Product Description        |
|------------|----------------------------|
| RTL2017CTN | 11.2" W x 10.2" L x 4.9" H |

- Insertable gel packs target up to 3 locations providing hot or cold therapy for relief of muscle inflammation and tension
- Ergonomically designed memory foam provides support for your head and neck while at home or on the go
- Soft velour cover enhances relaxation and is easily removable for cleaning
- Includes 3 re-usable gel packs and a travel bag for easy storage and portability

[Shop Now](#)



## Comfort Touch Gel Packs

**Item #**

RTL2017CTGP



[Shop Now](#)

## Comfort Touch™ Cooling Sensation Seat Cushion

**Item #**

RTL2017CTS

**Product Description**

18" W x 16" L x 3" H

- Ergonomic molded foam provides support and eases pain and discomfort associated with the pressure of prolonged sitting
- Gel top layer provides cooling sensation for added comfort
- Machine washable stretch cotton blend cover is low shear, water resistant, and comes with a non-skid bottom
- Molded foam cushion comes with a carrying handle to conveniently take cushion wherever you go




[Shop Now](#)



## Convolutated Foam Pad

| Item # | Product Description    |
|--------|------------------------|
| M6026  | 33" W x 74" L x 3.5" H |

- Convolutated foam mattress topper helps redistribute pressure and increase circulation
- Offers a cool, dry and comfortable surface for a more restful sleep

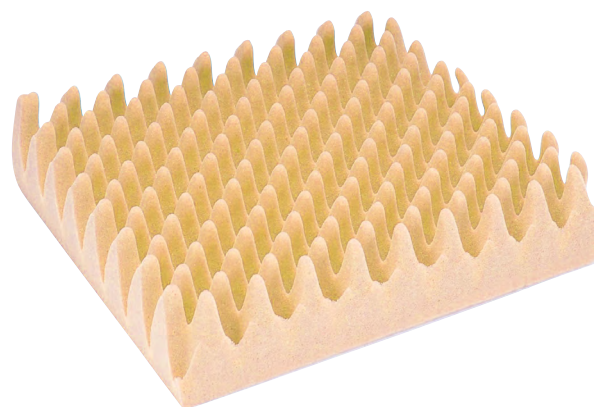[Shop Now](#)

## PrimeMat 2.0 Impact Reduction Fall Mats

| Item # | Product Description            |
|--------|--------------------------------|
| PM20BA | 24" W x 70" L x .75" H - Brown |
| PM20GA | 24" W x 70" L x .75" H - Grey  |

- Full length provides maximum coverage alongside bed
- Tough and lightweight - Mats are constructed from high quality molded Polyurethane foam inner core (210kgs/ cubic meter) to provide maximum absorption and rebound to reduce impact in the event of a fall
- Passed 5 foot "Egg Drop Test" - an egg was dropped from 5 feet above the mat without damage
- Low profile with tapered edges allow equipment such as overbed tables or patient lifts to roll over easily
- Anti-slip bottom netting keeps the mat secure to help prevent accidental tripping
- Embossed non-slip PVC leather look top is not only attractive but extremely durable and puncture resistant
- Non-slip pattern helps minimize slipping and is easy to clean with mild soap and water

[Shop Now](#)

## Safetycare Floor Mats 1 Piece with Masongard Cover

| Item #  | Product Description  |
|---------|----------------------|
| 7094    | 36" W x 66" L x 2" H |
| 7095-BF | 24" W x 66" L x 2" H |

- Safetycare floor mats help prevent injury in the event of an accidental roll-off from bed
- Made from high-density polyurethane foam



[Shop Now](#)

## Tri-Fold Bedside Fall Mat

| Item # | Product Description  |
|--------|----------------------|
| 14700  | 30" W x 72" L x 2" H |

- Helps to reduce the possibility of injuries from bed falls
- Conveniently folds in 3 sections for storage
- Made of high density foam to absorb impact from bed falls
- Non-skid bottom reduces the chance of the mat slipping
- Comes with durable, vinyl cover that is easy to clean



[Shop Now](#)



CORPORATE OFFICE 🏠 99 Seaview Blvd., Port Washington, NY 11050 📞 877-374-8763 🌐 [www.drivemedical.com](http://www.drivemedical.com)  
CANADA OFFICE 🏠 256 Aviva Park Drive, Vaughan, Ontario Canada L4L 9C7 📞 866-752-1363

©2025 Medical Depot, Inc. dba Drive DeVilbiss Healthcare, 99 Seaview Boulevard, Port Washington, NY 11050. All trademarks used in association with the sale of products of Drive DeVilbiss are trademarks owned by Medical Depot, Inc. All other trademarks, trade names, service marks and logos referenced herein belong to their respective companies. Refer to product instruction guide for warranty details.